


Facet Fixx®


Translaminar (TL)

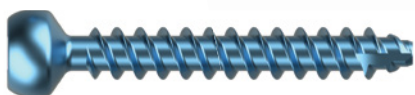

Cannulated
Ø4.5, 5.0mm
Lengths 35-60mm
(5mm Increments)

Lag feature maximizes compression
of the facet articulation process.

Transfacet (TF)

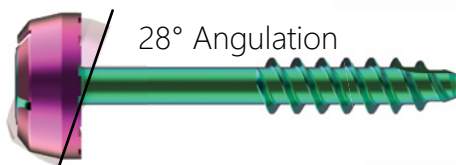


Cannulated
Ø4.5, 5.0mm
Lengths 25-50mm
(5mm Increments)

Washer undersurface designed to
prevent rotation without compromising
integrity of cortical surface.




Nexxt Spine, LLC
14425 Bergen Blvd, Suite B
Noblesville, IN 46060
(317) 436-7801
Info@NexxtSpine.com

70-021, Rev E

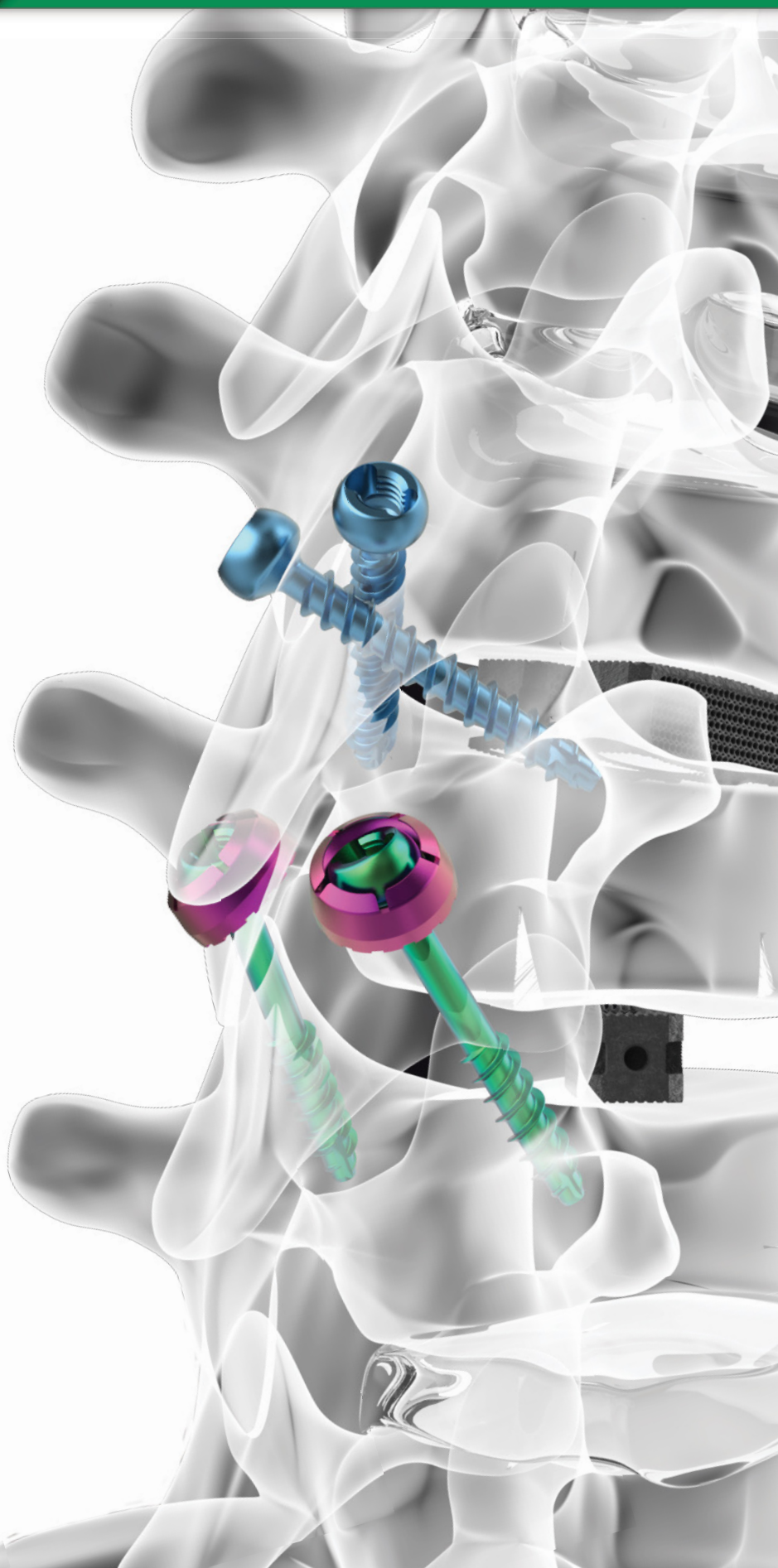

DTS Guide

Drill, Tap, Screw

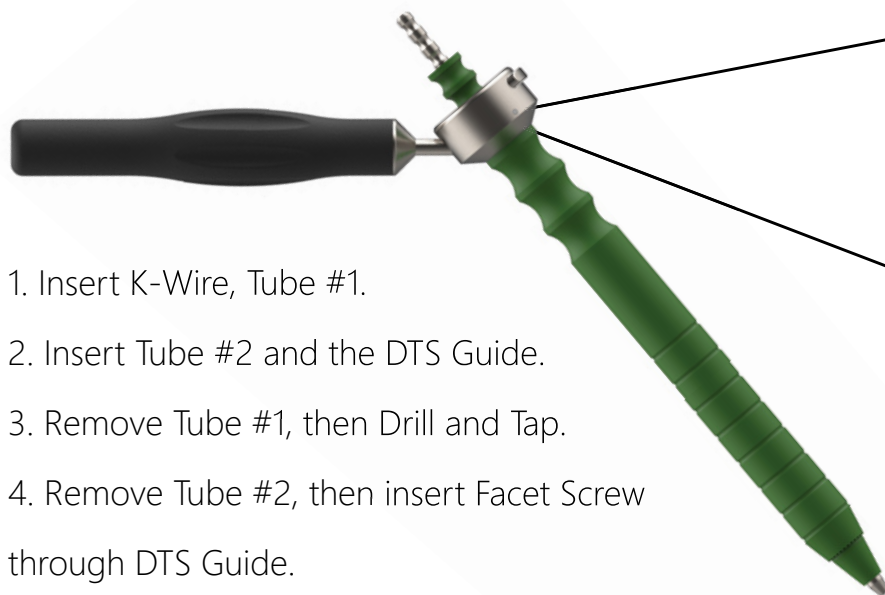

1. Insert K-Wire, Tube #1.
2. Insert Tube #2 and the DTS Guide.
3. Remove Tube #1, then Drill and Tap.
4. Remove Tube #2, then insert Facet Screw through DTS Guide.

Locking Feature retains
dilator tube during
procedure.

COMPATIBLE NEXXT MATRIX® PRODUCTS

Corpectomy



SA ALIF



ALIF




TLIF / rTLIF



TLIF Oblique



Lateral

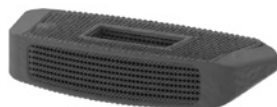


SA Cervical



Cervical



VBRs

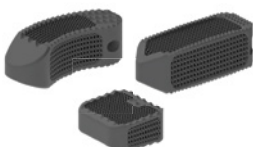

Nexxt Spine, LLC
14425 Bergen Blvd, Suite B
Noblesville, IN 46060
(317) 436-7801
Info@NexxtSpine.com
NexxtSpine.com



For indications, contraindications, warnings, precautions, potential adverse effects and patient counseling information, see the package insert or contact your local representative; visit NexxtSpine.com for additional product information.

All rights reserved. All content herein is protected by copyright, trademarks and other intellectual property rights owned by Nexxt Spine, LLC and must not be redistributed, duplicated or disclosed, in whole or in part, without the expressed written consent of Nexxt Spine, LLC. This material is intended for healthcare professionals, the Nexxt Spine sales force, and authorized representatives. Distribution to any other recipient is prohibited.