OVERWATCH

FEATURES AND BENEFITS

THORACOLUMBAR **FIXATION**

Amendia is now




Introduction/Implant Features

Power with Precision

The **Overwatch™ Thoracolumbar Fixation System** was designed to powerfully correct and stabilize the posterior thoracolumbosacral spine through its versatile collection of instruments and comprehensive implant offering.

The 5.5mm titanium rod based system offers a multitude of multi-axial pedicle screws, fixed and adjustable crosslinks and a variety of rod connectors and hooks. Multi-axial pedicle screws are available with standard or reduction tulip heads, self-tapping dual-lead and DualFix™ thread options, cannulated or non-cannulated and range from 4.5mm to 10.5mm in diameter.

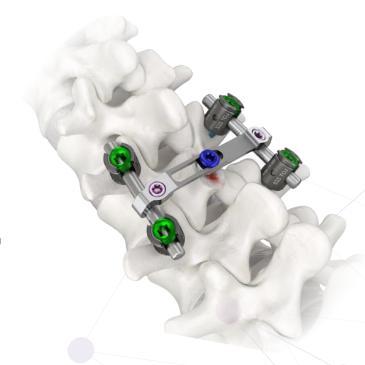



- Tulip groove feature for secure reduction tool engagement
- 60° of multi-axial angulation
- Constant tulip friction to maintain optimal head position
- Tulips with undercut buttress thread to minimize set screw cross-threading
- · Reduction Screws capable of 18mm of rod seating
- Available with dual-lead or DualFix™ threadforms

Power with Precision

Straight and Pre-Contoured 5.5mm Ti Alloy Rods



· Curved Pre-Cut Rods: 25mm - 250mm



· Straight Pre-Cut Rods: 30mm - 160mm



• Hex Rods: 300mm, 400mm, 500mm



Full assortment of hooks and rod connectors for construct adaptability





Adjustable and Fixed Crosslinks for increased construct torsional strength



Pedicle Screws, Set Screws, Cross-Links and Connectors designed with a hexalobe drive feature to minimize stripping during insertion

Instrument Features

Power with Precision

The **Overwatch™** Instruments were developed to be intuitive, precise and capable of performing in the most demanding scenarios. The system includes a full arsenal of refined tools for pedicle preparation, screw placement, rod contouring, reduction, rotation, compression and distraction.

